

Eliminate[™]

Aspiration catheter

THE DELIVERABLE SOLUTION
FOR THROMBUS REMOVAL

Eliminate™


Aspiration catheter

Eliminate™ is an innovative aspiration catheter designed to offer an optimal balance between crossing performance, kink resistance and thrombo-aspiration capability.

DELIVERABILITY


Eliminate's innovative technology provides improved deliverability and kink resistance.

- Pre-loaded stylet provides improved pushability and ensures optimal crossability in difficult to cross scenarios.
- Fully braided shaft construction for improved kink resistance of the catheter
- 23 cm long rapid exchange segment provides improved trackability.
- Hydrophilic coating facilitates easy navigation through tortuous anatomies.

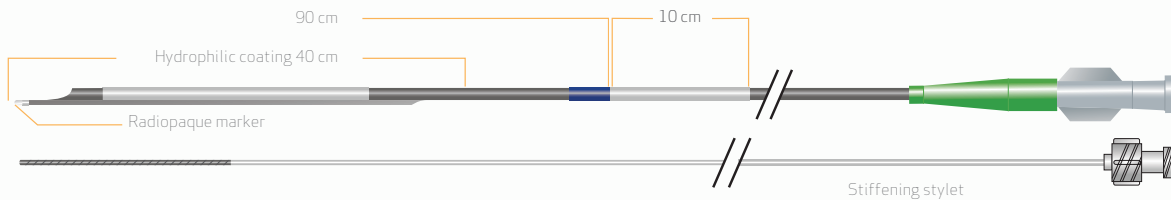

6 Fr guide catheter compatibility

Distal Part

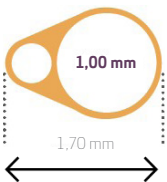

Extraction area
0,79 mm²


Proximal Part



Extraction area
0,95 mm²

7 Fr guide catheter compatibility

Distal Part

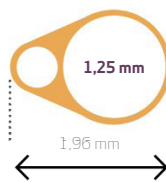

Extraction area
1,23 mm²

Proximal Part



Extraction area
1,33 mm²

THE DELIVERABLE SOLUTION FOR THROMBUS REMOVAL

ELIMINATE COMPONENTS



x2

Locking syringes
30 ml polycarbonate



x1

Extension line with
2-way stopcock



x1

Flushing tool



x1

Filter baskets
70 µm

PRODUCT SPECIFICATIONS

Tip Length:

6 mm

Coating:

Hydrophilic

Rapid exchange segment:

23 cm

Radiopaque marker band:

1 mm located at 4 mm from the tip

Marcas de posicionamiento:

10 cm = single white mark located at 90 cm
from distal tip

Wire compatibility:

0,014" (0,36 mm)

Catheter length:

140 cm

EASY TO USE

Single-handed operator

- Short and thin distal radiopaque marker provides excellent distal visualization and increases distal flexibility

• Complete kit contains:

- 30 ml locking syringes
- Extension line (25 cm) with 2 way stopcock
- Filter baskets

Fully braided
shaft



Stiffening stylet
with floppy distal segment



̳ 0,40 mm = 0,0157"

Tapered core wire

̳ 0,42 mm = 0,0165"

Floppy distal segment

Stainless Shaft

Eliminate™

Aspiration catheter

ORDERING INFORMATION

Eliminate	
Product Code	Guide catheter compatibility
EG1602	6Fr (D.I. - 0,70 pol / 1,78 mm)
EG1652	7Fr (D.I. - 0,80 pol / 2,03 mm)

LEARN MORE



**TERUMO
AMERICAS HOLDING, INC.**
2101 Cottontail Lane, Somerset,
New Jersey, 08873 U.S.A.
TEL: +1-732-302-4900
Tollfree (USA): +1-800-283-7866

**TERUMO MEDICAL
DE MEXICO S.A. DE C.V.**
Av. Insurgentes Sur 1647,
Piso 11 Col. San José Insurgentes
Del. Benito Juárez, C.P. 03900
México, D.F., México
TEL: +52-55-1085-0770

**TERUMO
COLOMBIA ANDINA S.A.S.**
Calle 100 #23-10 Piso 2
Bogotá D.C. Colombia Tel:
+57-1256-8400

**TERUMO
LATINAMERICA CORPORATION**
Doral Corporate Center I,
8750 NW 36th Street, Suite 240,
Miami, FL 33178 U.S.A.
Tel: +1-305-477-4822

**TERUMO
MEDICAL DO BRASIL LTDA.**
Praça General Gentil Falcão, No.
108-9º andar
Cidade Monções - Brooklin Novo
CEP: 04571-150
São Paulo - SP, Brasil
TEL: +55-11-3594-3800

TERUMO CHILE LTDA.
Avenida Apoquindo 2929
Piso 16 Las Condes Santiago
de Chile 7550246
TEL: +56-2-480-9600



©TERUMO CORPORATION. Jun, 2021.
All names are trademarks or registered trademarks of TERUMO
CORPORATION and its proprietary owners.

TERUMO CORPORATION
2-44-1 Hatagaya, Shibuya-Ku, Tokyo, 151-0072 Japan www.terumo.co.jp